

QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)

TURNING ON/OFF, CHARGING THE DEVICE




On/off
button on
right side



Power
supply
connector

START MENU

Note: this screen (test options) can vary depending on your set up

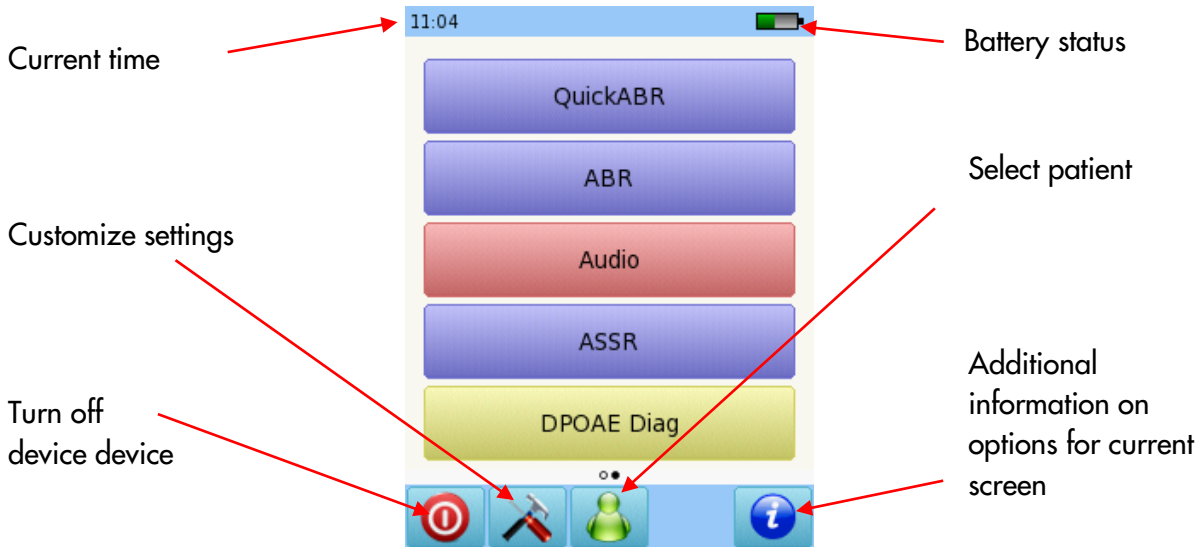

QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)

TESTING

AUDITORY STEADY STATE RESPONSE (ASSR)

Connect the electrode cable



to the white connector.




Connect the transducer to the appropriate connector:

Headphones or insert earphones: blue connector

Probe: red connector

Second probe for binaural testing: blue connector

Clean the patient's skin and place the electrode pads. Connect the clips to the corresponding sites.
Recommended montage is:

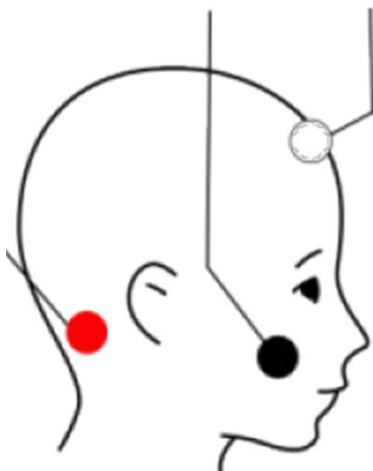

White: high forehead

Black: cheek

Red: mastoid

For binaural testing, place red electrode on the nape.


QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)

If you are using an ear probe or insert earphones, select the appropriate ear tip size and place it on the probe tip. Insert the probe in the ear making sure that the ear tip fits tightly into the ear canal.


If you are using headphones, place them over the ears.

Select ASSR

Select the protocol you would like to use. To change the parameters, select one of the presets and then the tools button at the bottom of the page:

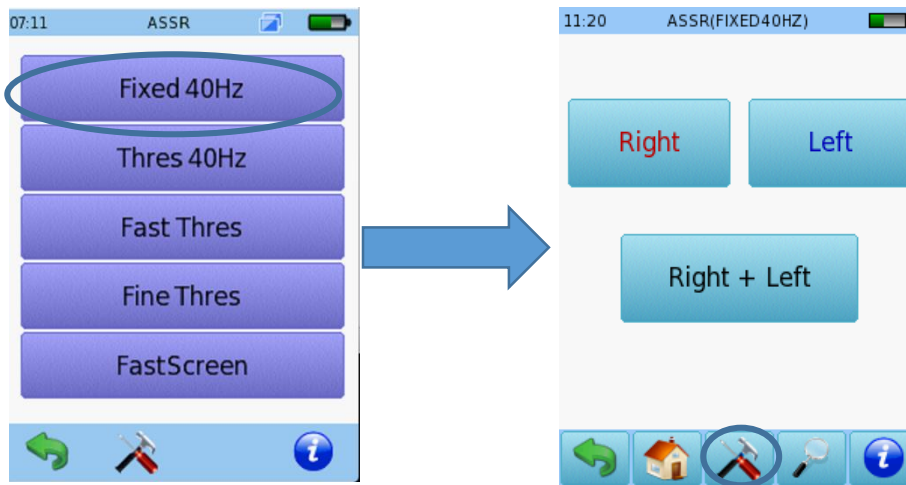

CONFIGURABLE TEST PARAMETERS

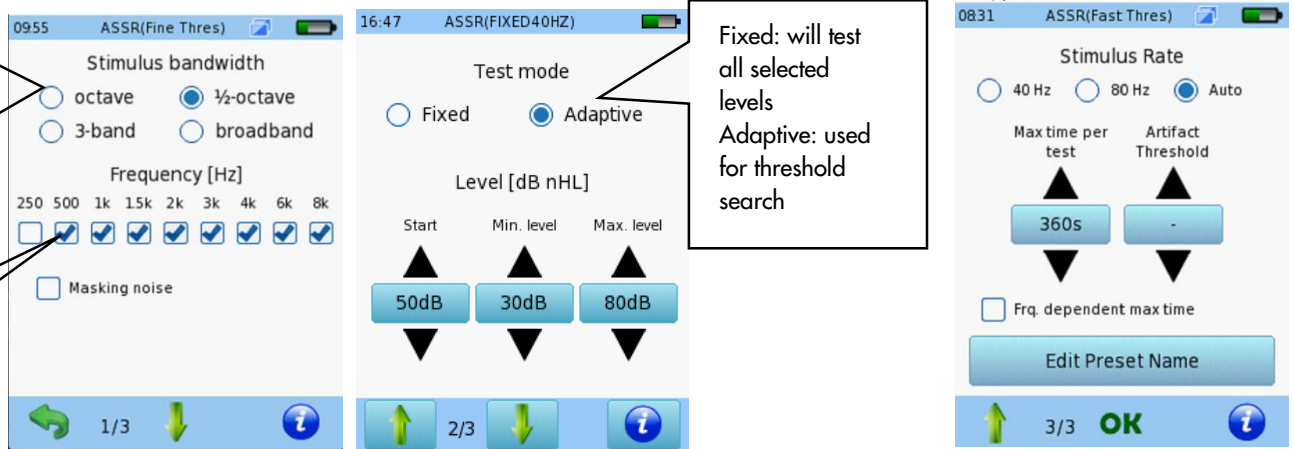

Half octave stimuli provide better frequency specificity but longer test time

Select frequencies to test

Select masking noise (monaural testing only)

Select start, minimum and maximum level for adaptive mode

Fixed: will test all selected levels
Adaptive: used for threshold search

40 Hz: awake and adult patients, low frequencies ≤ 2 kHz
80 Hz: babies and higher frequencies
Auto: mixed rate, optimized for performance


QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)

Select the test ear by pressing either the Right or Left button, or Right + Left for binaural testing. Binaural testing is available when headphones, insert phones or two ear probes are connected.




Make sure the patient is sitting or lying down comfortably.

An impedance check will be performed prior to measurement. If impedance is below 6 k ohm, the start button will become available. You can configure the settings for the test to *autostart* if the impedance is lower than 4 k ohm. If impedance is higher than 6 k ohm but below 12 k ohm, a skip button will become available. An impedance below 6 k ohm is recommended for best test results.

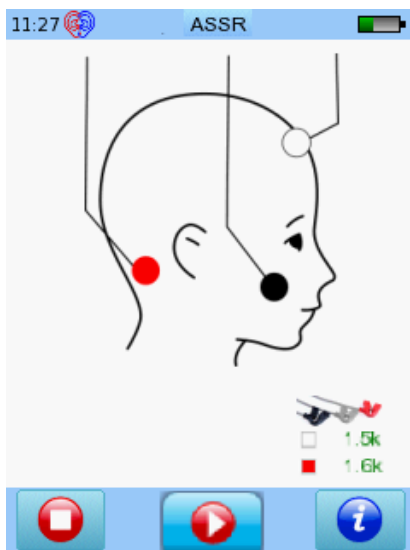


QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)

If you are using an ear probe, a calibration will be performed together with the impedance check. You can select the simple calibration or the expert calibration mode in the device set up.

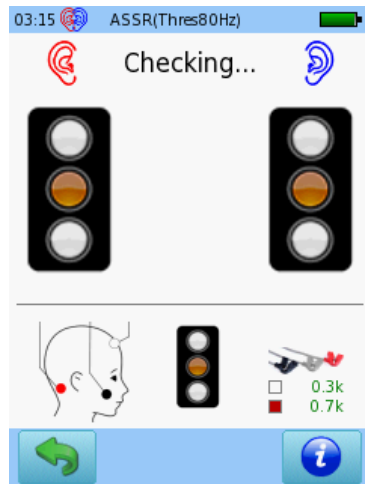

Simple calibration mode

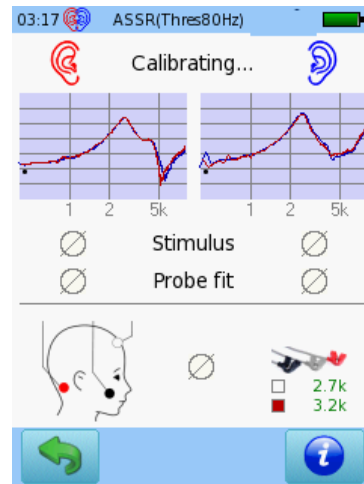

Expert calibration mode

ASSR TEST SCREEN

Annotations for the ASSR TEST SCREEN:


- Progress for each tested frequency
- Frequency being tested
- Intensity level being tested
- Maximum test time
- Overall EEG value should remain in the lower third. Press to enable raw EEG view
- Impedance monitor
- Progress
- Stop test
- Pause ongoing test
- Skip this single test

QUICK START GUIDE SENTIERO ADVANCED ASSR (GOLD LABEL)


ASSR TEST RESULTS

Please note that the estimated audiogram is only available in the adaptive test mode.



Note: These circles indicate that other screens are available to view. Swipe your finger right or left to get other views.

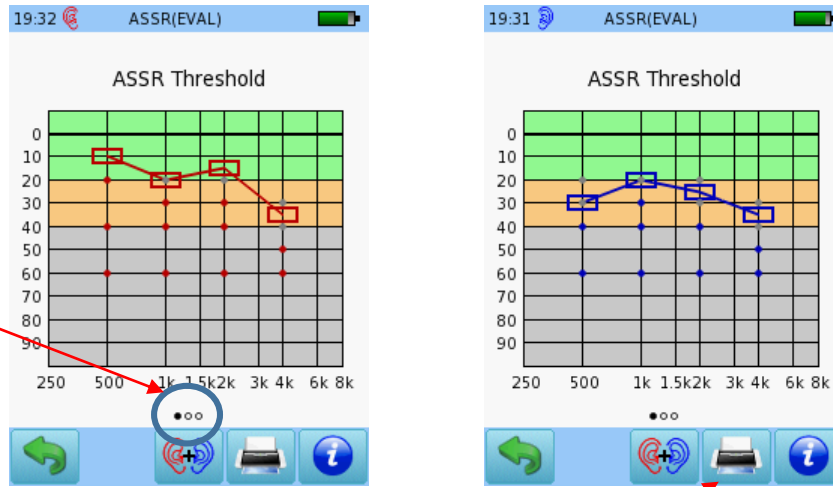

Print results to a label or to a pdf file

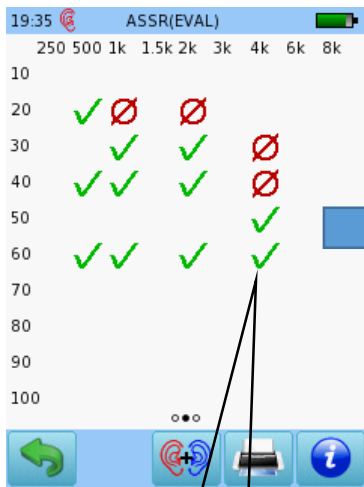

Press on a check mark or circle to view details

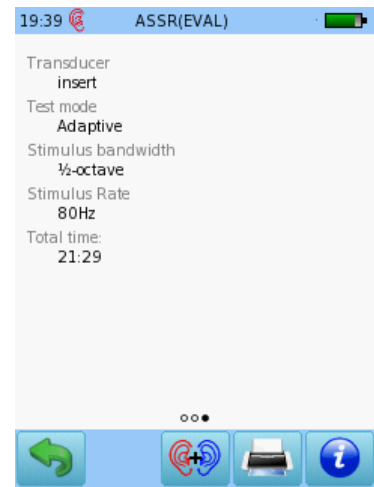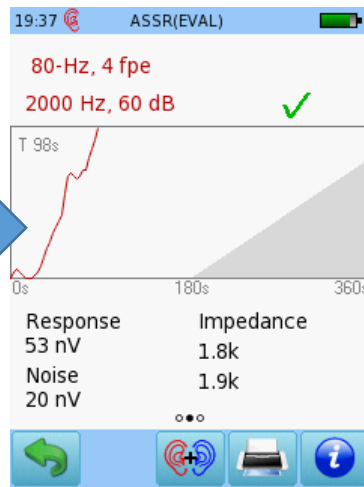

QUICK START GUIDE

SENTIERO ADVANCED ASSR

(GOLD LABEL)

Quick Guide Information

Release date: 06/2022
Revision: 02
Valid from: Firmware Rev. 2.8, Mira PC Software Rev. 2.4
Applies to: Model: SOH100360

All mentioned items, products, brands and trademarks are registered or owned by the mentioned companies.

All information, illustrations, and specifications provided within this manual are based on the latest product information available at the time of publication. PATH MEDICAL reserves the right to make changes at any time without notice.

Please note that not all products displayed on our webpage or advertising materials are available in all territories. Please contact your local distributor or PATH MEDICAL for further information.

The latest revision of the Quick Start Guide is available online at www.pathme.de/download.

Errors and omissions excepted.

Copyright Notice

No part of this manual may be reproduced, translated, stored, or transmitted, in any form or by any means, electronic, mechanical, photocopying, recording, or otherwise, without the prior written consent of PATH MEDICAL GmbH.

Copyright © 2022 PATH MEDICAL GmbH

PATH MEDICAL GmbH
Landsberger Straße 65
82110 Germering
Germany

Tel.: +49 89 800 76 502 Fax: +49 89 800 76 503 Internet: www.pathme.de

CE
0124

